




Disposable Sampling Swab

Product Description

96000J is a Flocked Sampling Swab used for cell & virus specimen collection, and it has been designed for Nasopharyngeal clinical diagnostics. It utilizes state of the art “spray on technology” that the flocking process by means of an electro-static charge perpendicularly attaches millions of nylon microfibers on the medical grade handle tip. The flocked swab is ideal for collecting large amount of cells and rapid elution of the specimens that instantly releases the cells into the transport medium. It has been well recognized and adopted by the diagnostic test kit manufacturers who produce reagents in molecular genetics, forensics and clinical laboratories sectors.

Product Features

- Ergonomic and anatomic design
- Improved sample collection
- Superior sample elution
- Increased assay sensitivity
- Quantitative volume transfer
- Certified free of inhibitors and interference

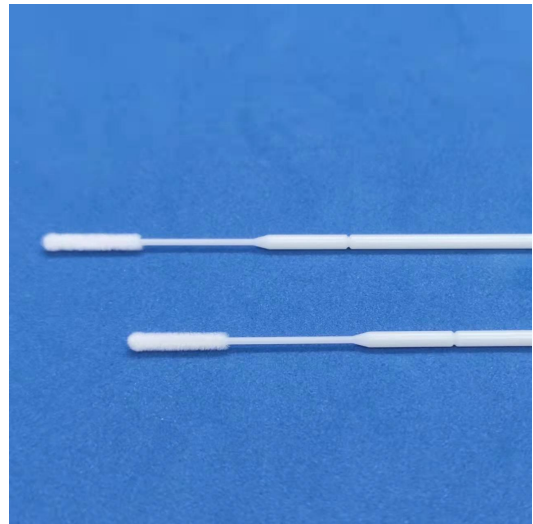

Product Specifications

| Item No. | Dimension of flocked tip | | | Dimension of ABS handle | | | |
|--------------------------|--------------------------|-----------|--------|--|--------------|------------|--------------|
| | Width | Thickness | Length | Diameter 1 | Diameter 2 | Breakpoint | Total length |
| 96000J | 3±0.5mm | 3±0.5mm | 16±3mm | 2.5 mm | Taper to 2.0 | 50±3mm | 148±15mm |
| Packaging | | | | 100 pcs/bag, 50 bags/ctn, 5000 pcs/ctn | | | |
| Sterilization method | | | | Ethylene oxide sterilization | | | |
| Ethylene oxide residue | | | | ≤10µg/g | | | |
| Best storage conditions | | | | 4°C— 30°C | | | |
| Quality guarantee period | | | | 3 years | | | |



Disposable Sampling Swab

Directions

Step 1: Put the tip of swab into nose.

Step 2: Squeeze with your nose for about 1 minute to allow the swab to completely absorb sampling.

Step 3: Put the swab on test machine.



Miracle Technology Co.,Ltd.

Address: Room 301, Building A, No.18, Rongshuxia Industrial Zone, Tongxin Community, Baolong Street, Longgang District, Shenzhen, China, 518116

Tel: +86-0755- 89616773 Fax: +86-755-89616775

Web: www.mantacc.com



Share Info Consultant Service LLC Repräsentanzbüro

Address: Heerdter Lohweg 83, 40549 Düsseldorf

Tel: 0049 1767 0057022

Fax: 0049 1767 0057022

Emai: info@share-info.cn

DIMDI NO.: DE/0000047946

CE Logo

