




Disposable Sampling Swab

Product Description

93050M is a Foam Sampling Swab used for cell & virus specimen collection, and it has been designed for oropharyngeal clinical diagnostics. This specialty applicator is tipped with 100ppi medical grade polyurethane foam. It has been well recognized and adopted by the diagnostic test kit manufacturers who produce reagents in molecular genetics, forensics and clinical laboratories sectors.

Product Features

- Medical grade quality
- Excellent collection and release properties
- Soft foam tip for patient comfort
- Used dry or with media
- Made in China

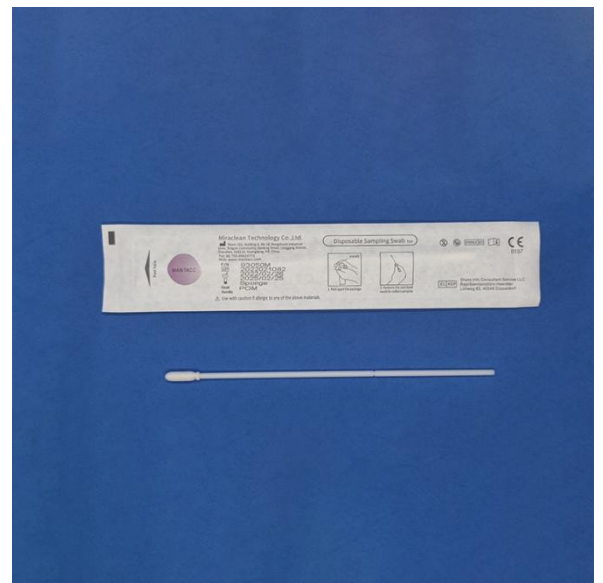

Product Specifications

Item No.	Dimension of foam tip			Dimension of POM handle			
	Width	Thickness	Length	Diameter <u>1</u>	Diameter <u>2</u>	Breakpoint	Total length
93050M	4.4±0.5mm	4.4±0.5mm	15±3mm	2.5 mm	Taper to 2.0	95±3mm	152±15mm
Packaging				100pcs/bag, 50bags/ctn, 5000pcs/ctn			
Sterilization method				Ethylene oxide sterilization			
Ethylene oxide residue				≤10µg/g			
Best storage conditions				4℃— 30℃			
Quality guarantee period				3 years			



Disposable Sampling Swab

Directions

Step 1: Put the tip of swab into mouth.

Step 2: Squeeze with your tongue for about 1 minute to allow the swab to completely absorb saliva.

Step 3: Put the swab on test machine.



Miraclean Technology Co.,Ltd.

Address: Room 301, Building A, No.18, Rongshuxia Industrial Zone, Tongxin Community, Baolong Street, Longgang District, Shenzhen, China, 518116

Tel: +86-0755- 89616773 Fax: +86-755-89616775

Web: www.mantacc.com



Share Info Consultant Service LLC Repräsentanzbüro

Address: Heerdter Lohweg 83, 40549 Düsseldorf

Tel: 0049 1767 0057022

Fax: 0049 1767 0057022

Emai: info@share-info.cn

DIMDI NO.: DE/0000047946

CE Logo

